



MONITORING



PROTECT CRITICAL EQUIPMENT



COMPLIANCE



NOTIFICATION




GRAPH OF EVENTS

Continuous Temperature Monitoring

New RTD Temperature Sensor

The same company that has provided protection and peace of mind for over 30 years is introducing their 3 wire precise and accurate Resistance Temperature Detector sensor. Our new RTD sensor comes equipped with a high accuracy wide range probe and one sensor covers the full laboratory temperature range.

- * Equipped with a PT100 platinum RTD probe
- * Operates over a -200° C to +100° C range
- * Offers a .07 degree resolution
- * Single or 2 input models
- * Standard version has typical 0.4 degree accuracy over span. Certified version has 0.1 degree accuracy over span (powered).
- * Certified version comes with a certificate showing 0.1 degree accuracy at -196, -70, -40, -20, 0, 20, and 50 degrees C.
- * RTDs offer greater stability and accuracy than thermocouples or thermistors.
- * Sensor operates powered or with battery

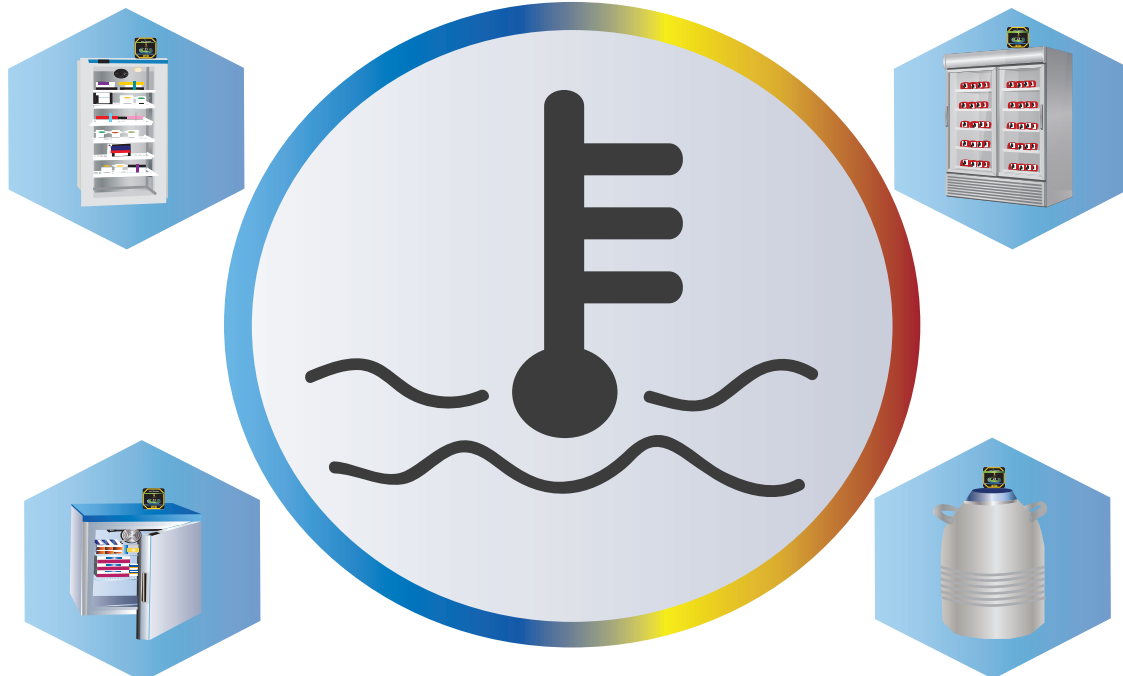



MONITORING



PROTECT CRITICAL EQUIPMENT



COMPLIANCE



NOTIFICATION



GRAPH OF EVENTS

Continuous Temperature Monitoring

How it works

