


Rees Scientific Dialer Options



Back View

CELLULAR

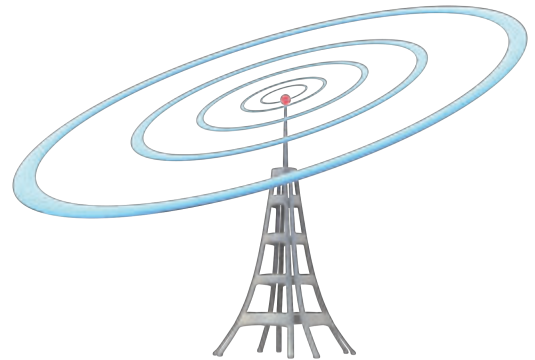

Rees Scientific has improved and expanded their capabilities and offer more options for their alarm dialing.

VOIP Dialer

- * For use with virtualized Centron Presidio installation
- * Constant connection to box
- * Less reliant on traditional analog solutions
- * Highest level of security for mission critical assets and events
- * Optional 1 or 4 line dialer

VOIP Software

- * Viable option for sites that don't own or want POTS line and TAPI modem
- * Ability to leverage modern network infrastructure
- * Less hardware required (eliminates need for modem, telephone line and internal connection of modem)

VOIP Backup Dialer

- * Monitors phone line for connectivity
- * Additionally, if the main phone line service is interrupted, the VOIP backup dialer will take over alarm duty
- * Another safeguard ensuring your critical equipment is safe and protected

Cellular Dialer

- * External Cellular Dialer to connect to customer provided PC
- * Performs telephone alarm dial-out for the Rees Scientific Virtual Node and requires 1 external USB port
- * Synthesized voice dial out capable of working with RSC customizable voice tags
- * Comes with prepaid GSM SIM Card, but ongoing cellular service through a Cellular GSM telephone provider is required
- * Can be main dialer or back up