


Feel safe and secure through the Rees Scientific Automatic Failover Node Protection

You now have the ability to protect your entire Centron Monitoring System against hardware failure through the Automatic Failover Node. The Automatic Failover Node continuously monitors your main node and takes over if a problem occurs. This protects against failure of hard drives, power supplies, phone lines, dialers, network connection failure, etc. It also protects against failure of inter-site communications. If a remote site or sites fall out of contact with the main site, the Automatic Failover Node(s) at the remote site(s) will take over and provide full data logging and alarm handling functions.

- Monitor up to 4 nodes with one machine
- Automatic Failover Node(s) can be located in the Server Room or in any workstation
- Main node can be any Centron node type

Ideal for:

- * All critical systems where service interruption could cause serious issues
- * Systems with one or more remote sites

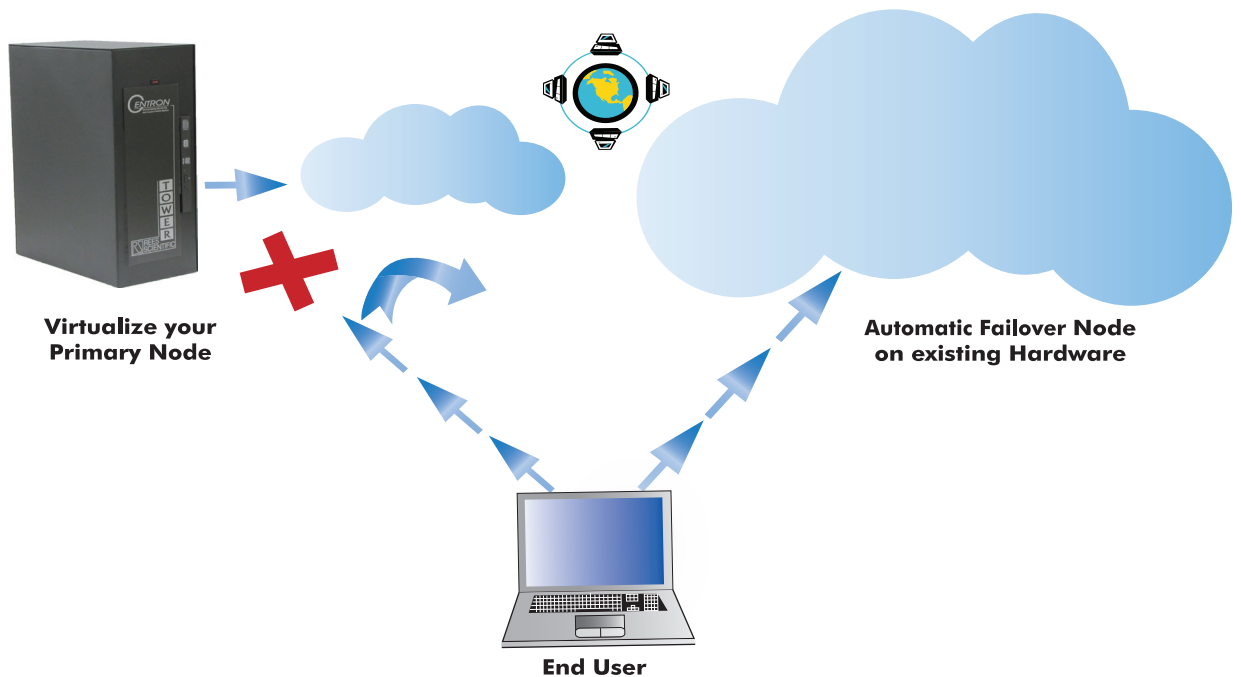

PRNT2528